

AUTOMATED INVENTORY MANAGEMENT COMPACT INTERNATIONAL RFID REFRIGERATOR - TS060



COMPACT RFID REFRIGERATOR - TS060

RELIABLE INVENTORY SOLUTIONS:

- + Increase product visibility and reduce shrinkage
- + Eliminate inefficient manual reconciliation processes

IMPROVE REGULATORY COMPLIANCE:

- + Generate policy compliance
- + Remote temperature monitoring

ABSOLUTE ACCOUNTABILITY:

- + Ensure product integrity
- + Proactive alarms for expiring products
- + Integrate with existing inventory applications

INVENTORY OPTIMIZATION:

- + Usage tracking for improved charge capture rates
- + Prevent stock outs

PERFECT FOR USE IN:

- + Biobanks
- + Tissue storage facilities
- + Hospitals
- + Research labs
- + Pharmacies
- + Stockrooms

SPECIFIC APPLICATIONS INCLUDE:

- + Reagents
- + Bio-samples
- + Pharmaceuticals
- + Cardiovascular products
- + Orthopedic supplies



AN INTERNET-ENABLED RAIN RFID ENCLOSURE IS PLACED AT THE POINT OF USE IN YOUR FACILITY.



PRE-AUTHORIZED USERS ACCESS THE ENCLOSURE AS NEEDED.



WHEN A USER ADDS OR REMOVES AN ITEM, THE SYSTEM AUTOMATICALLY TRACKS THE TRANSACTION IN REAL-TIME.



ALL TRANSACTIONS ARE RECORDED IN REAL-TIME AND ACCESSIBLE VIA JETSTREAM, TERSO'S CLOUD API PLATFORM.

Learn how Terso Solutions can free you from the complications and inefficiencies of traditional inventory management systems. Call +49.621.8501.233 or visit WWW.TERSOSOLUTIONS.COM

Traditional inventory management systems are slow, complex and just plain frustrating to work with. On top of that, an outdated system takes precious time away from doing what you do best.

TERS0 SOLUTIONS CAN FREE YOU FROM ALL THE CHAOS AND CONFUSION.

Our RAIN RFID-enabled enclosures automatically track and monitor inventory 24/7 which eliminates manual processes, increases visibility, and captures information in real-time. The inventory transactions can be integrated into existing enterprise resource planning (ERP) and inventory management workflows.



TERS0TM
SOLUTIONS, INC.



**LIBERATE YOURSELF FROM
THE COMPLICATIONS OF
AN OUTDATED INVENTORY
MANAGEMENT SYSTEM.**

**COMPACT RFID REFRIGERATOR - TS060
104 L INTELLIGENT REFRIGERATOR**

PHYSICAL CHARACTERISTICS

CHARACTERISTIC	Metric (mm) H x W x D
External Dimensions	1010 x 595 x 620
Internal Dimensions	618 x 492 x 415
External Space	1010 x 595 x 1270, Place 76.2mm From Back Wall
Volume	104 L
Weight	75 kg
Security Device	Electromagnetic
Refrigerant	R - 134 A
Cooling Method	Forced-air
Temperature Range	2°C to 12°C

ENVIRONMENTAL


Environmental	IP41, RoHS Compliant
---------------	----------------------

COMMUNICATIONS

RFID Frequencies	864 to 868 MHz
Communications	Ethernet, Cellular, Optional WiFi
Access Control	LF Proximity Reader

POWER

Electrical Power	230 VAC, 50 Hz, 2A and 106W
Power Consumption	84W max
Battery	Two - 12VDC, 1.3Ah, Sealed Lead Acid

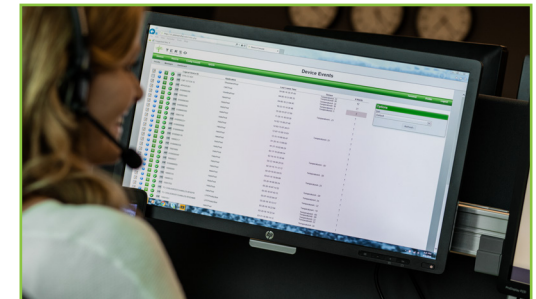

REFRIGERATOR FEATURES

- Open shelves
- Configurable drawers
- Remote temperature monitoring
- Custom branding
- Ergonomic riser

